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### Instructions for use

## Friction–Fit–System

Before using this product, please read the instructions for careful use.

The manufacturer will not accept any liability for damage resulting from non–compliance with these instructions for use.

#### 1. Description of the product

The friction fit system (FGP) represents a reliable possibility to achieve permanent friction for a newly produced dental appliance or in case of lost metal fit. The friction acrylic allows to obtain uniform friction which can be individually adjusted depending on the degree of adhesion.

#### 2. Indication

Telespic and conical crowns, intra- and extra-coronal attachments but also implant restorations can be reactivated with this system. Additionally, FGP offers the possibility to produce secondary elements of telescopic and conical as well as CoCr supply according to the less expensive single-piece casting technique.

#### 3. Contraindication

No more than 2 crowns can be rebased for direct use in situ. If more crowns are rebased, it may no longer be possible to remove the denture or teeth may be extracted.

#### 4. Hazard and safety information

##### FGP friction acrylic

##### Signal Word:

Warning

##### Hazard statements:

H315	Causes skin irritations.
H317	May cause an allergic skin reaction.
H319	Causes serious eye irritation.
H335	May cause respiratory irritation.

##### Precautionary statements:

##### Prevention:

P273	Avoid release to the environment.
P280	Wear protective gloves/eye protection.

##### Response:

P302+P352	IF ON SKIN: Wash with plenty of soap and water.
P305+P351+P338	IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.

##### FGP–bonding agent

##### Signal Word:

Danger

##### Hazard statements:

H225	Highly flammable liquid and vapour.
H319	Causes serious eye irritation.
H336	May cause drowsiness or dizziness.
EUH066	Repeated exposure may cause skin dryness or cracking.

##### Precautionary statements:

##### Prevention:

P210	Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking. Keep container tightly closed.
P233	Avoid breathing vapours.
P261	Avoid breathing vapours.
P280	Wear protective gloves/protective clothing/eye protection/face protection.

##### Response:

P303+P361+P353	IF ON SKIN (or hair): Remove/Take off immediately all contaminated clothing. Rinse skin with water/shower.
P305+P351+P338	IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.

##### Storage:

P403+P235      Store in a well–ventilated place. Keep cool.

**Product contains:** butanone

##### FGP insulation agent

##### Signal Word:

Danger

##### Hazard statements:

H225	Highly flammable liquid and vapour.
H319	Causes serious eye irritation.
H336	May cause drowsiness or dizziness.

**Product contains:** propan–2–ol

#### 5. Storage and durability

##### 5.1 Storage information

FGP bonding agent and FGP components A and B must be stored dry at a temperature of 4 to 25° C (39 to 63 °F) and protected against light.

##### 5.2 Durability information

The shelf life printed on the packaging is guaranteed if proper storage is ensured. Once the shelf life provided on the container has been expired, the products must not be used any longer.

#### 6. Proper disposal

Please adhere to the information provided in the safety data sheets of the products.

#### 7. Cleaning

For cleaning the denture it is recommended to use toothpaste and a toothbrush and no commercially available cleaning tabs.

#### 8. Processing

##### 8.1 New production of a telescopic or conical restoration in the single–piece casting technique

Produce the primary elements as usual. Produce a coping from thermo–forming foil or immersion wax with a thickness of at least 0.2 mm which is used as a spacer for the FGP acrylic. Cut off this coping approx. 1 mm above the cervical margin so that the cervical margin of the secondary element can be finished thinly later on. Then duplicate the model with the primary elements. The chrome cobalt structure and the secondary elements are modelled onto the investment material model. After casting, finish and polish the CoCr frame. If the secondary elements are to be veneered with heat–polymerizing resin, this work needs to be carried out prior to the FGP relining. The FGP acrylic will not resist these temperatures. Then the actual FGP relining work can be started.

##### 8.1.1 Preparation of the object

Blast the surface onto which FGP acrylic has been applied. Please use blasting nozzles for manual sandblasting. Blast with 50 µ aluminum oxide at 2 to 3 bar (units with automatic blasting material circulation are not suitable due to the contamination of the blasting material; they contaminate the metal surface so that the adhesion of the FGP bonding agent and the acrylic is no longer ensured).

##### Please note:

Blasting material residues are removed from the blasted metal surface using oil-free compressed air or a clean, dry brush. Uniform application of FGP bonding agent is required to achieve proper bonding. Avoid the formation of puddles and visible brush strokes. Then the FGP bonding agent is applied one time using a clean brush and dried for approx. 30 sec. The bonding agent is exclusively intended for the special FGP metal fit (non-precious metal and precious metal).

##### 8.1.2 Mixing

Mix the acrylic components A + B in the ratio 1 : 1 and stir the acrylic material without the formation of bubbles. The acrylic components A + B harden 180 sec after the beginning of the stirring process. Pour the material into the secondary elements; pour slightly more material into the occlusal area so that the acrylic will flow without forming any bubbles. In order to control the friction, FGP insulating agent can be applied onto the primary element. The more insulating agent is applied, the lower the friction will be. The acrylic hardens after approx. 420 sec. Any insulating residues must be removed from the denture by holding it under running water for a short time. The insulating agent is water-soluble.

##### 8.2 Repair of a worn-out metal-to-metal fit

First grind the surface to be relined with an aluminium-oxide-bonded stone or a new diamond so that a space of at least 0.2 mm is obtained. Then continue as described under item 8.1.1. Before starting with item 8.1.2, first the dentist has to remove tartar and saliva from the primary elements in the patient's mouth and apply liquid vaseline to insulate these parts. Prior to this, it is recommended to prepare a model of the oral situation and to check for undercuts.

##### Please note:

FGP insulating agent must not be used in the mouth. In such cases, vaseline must be applied to insulate the primary elements. With the help of the individual application of liquid vaseline the layer thickness can be increased (the thicker the layer, the lower the friction!). Use impression silicone to block out rough spots or undercuts on primary elements (also in the margin area). In cases of severe gingival retraction or swollen gingiva, a retraction cord must be placed. The object, which has been filled with FGP acrylic material, is placed onto the primary elements and pressed onto them manually. The patient should exert slight masticatory pressure for precise adjustment. Acrylic material that has flown over the cervical margin must be removed. Approx. 90 to 120 seconds after starting the mixing process the denture should be removed from the primary elements and placed on again. The FGP acrylic material has hardened after 7 minutes; the impression can now be removed and the margins can be trimmed.

**Please note:**
If undercuts at primary elements – which are covered by the outer telescope – are not noticed, it is no longer possible to remove the denture; therefore it is necessary to block out such spots! In the case of cover dentures, remove the denture acrylic in the vestibular area of the secondary elements to be rebased. Only this way it is possible to remove acrylic material that has flown over in the cervical area. If it is not removed and hardens, problems may result when removing the denture.

#### 8.3 Alternatively

The secondary elements can be rebased on an exact model made of super-hard stone, model resin or on milled metal dies. For this purpose an individual tray must be produced first. If the friction that is achieved in this way is not sufficient, rebasing must be repeated in the mouth using a minimum quantity of FGP acrylic. The components A + B are processed as described under item 8.1.2.

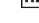
#### 9. Sources of errors and how to eliminate them

FGP friction acrylic hardens too quickly:	Temperature of the secondary elements too high during pou ring in the FGP acrylic or the FGP acrylic was too warm. Repeat entire work process.
---	--

Bubbles in the FGP friction acrylic after Simply mix a few drops of acrylic, pour in and the removal:	place onto the insulated die.
---	-------------------------------

#### 10. Symbols

##### 10.1 Symbols according to DIN EN ISO 15233–1

	Reference number
	Manufacturer
	Lot number
	Can be used until YYYY–MM–DD
	Temperature limits, upper limit: 25° C, lower limit: 4° C
	Dry storage
	Attention!

##### 10.2 Symbols according to GHS regulation

	GHS 02
	GHS 07

##### 10.3 Symbols according to NFPA (USA)

		Hazard diamond
---	---	----------------

#### 11. Additional information

These instructions for use are based on state-of-the-art methods and equipment and our own experience. The product may only be used for the indication described under item 2. The user himself is responsible for processing the product. Liability for incorrect results shall be excluded since the manufacturer does not have any influence on processing. Any occurring claims for damages may only be made up to the value of our products.

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